

Zimmer Biomet Dental Tapered Osteotomes

The Tapered Osteotomes are available in both straight and angled designs. There are four straight and four angled Tapered Osteotome sizes available. The sizes correspond with the diameters of Zimmer Biomet Dental Tapered Implants (3.25 mm, 4.0 mm, 5.0 mm and 6.0 mm diameters).

The depth markings are machine etched into the matte finished tip surface. This enhances depth line visibility, makes the depth line permanent and reduces the possibility of corrosion on the operating tip.


Angled Design Note: Angled design has a 30° Offset



Straight Design




Striking end of the Osteotome

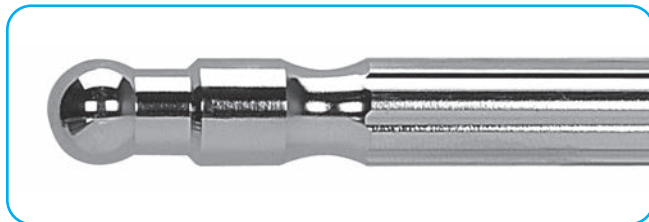

The ball at the striking end of the osteotome is forgiving of off-angled mallet strikes, delivering force in a linear vector to the operating end.

Operating end of the Osteotome

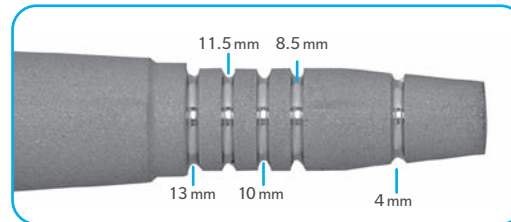

The depth lines are aligned to indicate crestal placement of the Certain and Ex Hex Tapered Implants.

Straight And Angled Osteotomes

Description	Item No.
Straight Tapered 3.25 mm	NTOST3
Straight Tapered 4.0 mm	NTOST4
Straight Tapered 5.0 mm	NTOST5
Straight Tapered 6.0 mm	NTOST6
Angled Tapered 3.25 mm	NTOST3A
Angled Tapered 4.0 mm	NTOST4A
Angled Tapered 5.0 mm	NTOST5A
Angled Tapered 6.0 mm	NTOST6A
Straight Tapered Kit (Includes all four sizes and sterilization case)	NTOST0
Angled Tapered Kit (Includes all four sizes and sterilization case)	NTOST0A



Standard Osteotome Procedures Are To Be Followed When Using The Tapered Osteotomes.

Unless otherwise indicated, as referenced herein, all trademarks are the property of Zimmer Biomet; and all products are manufactured by one or more of the dental subsidiaries of Zimmer Biomet Holdings, Inc. and marketed and distributed by Zimmer Biomet Dental and its authorized marketing partners. For additional product information, please refer to the individual product labeling or instructions for use. Product clearance and availability may be limited to certain countries/regions. This material is intended for clinicians only and does not comprise medical advice or recommendations. Distribution to any other recipient is prohibited. This material may not be copied or reprinted without the express written consent of Zimmer Biomet Dental. ZB1087 REV A 04/20 ©2020 Zimmer Biomet. All rights reserved.

